


### EMG Accessories

Our range of EMG electrodes has been extended to include various disposable and reusable electrodes. By keeping a wide selection of functional diagnostic electrodes for EMG studies we trust that you will find the appropriate electrodes for your specific study and patient.

Our electrodes are:

- Specifically designed for pelvic floor diagnostics
- Practical and innovative

- Disposable - eliminating the risk of cross-contamination
- Designed for enhanced patient comfort through practical solutions

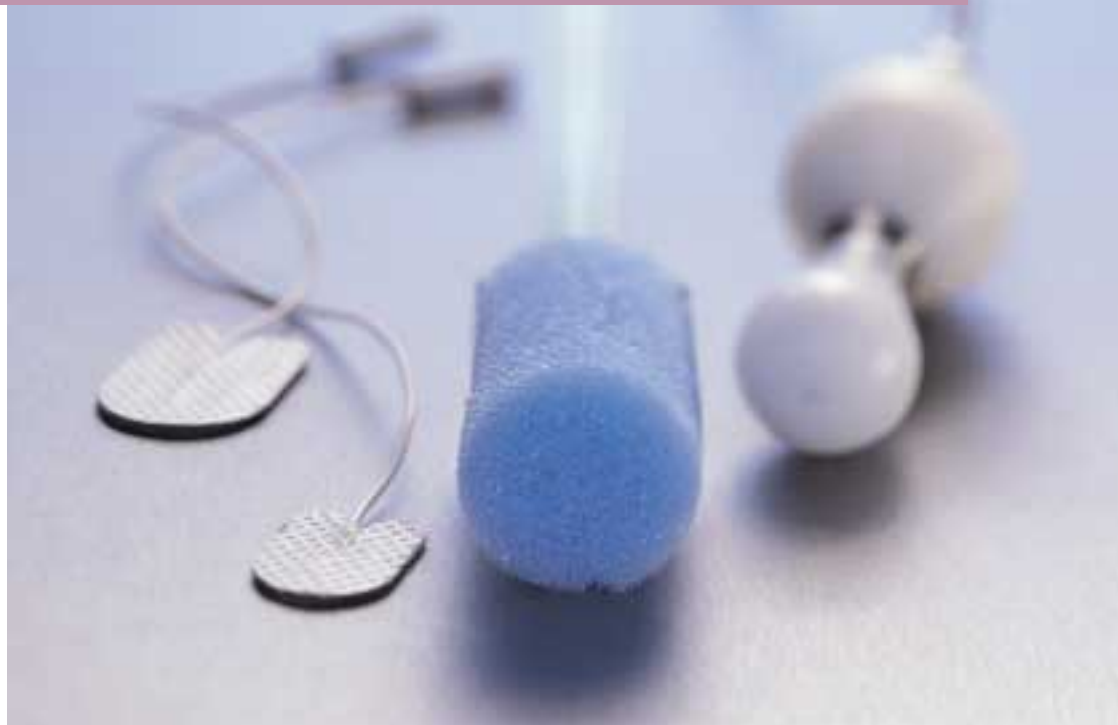


## EMG Accessories

### Surface Electrodes

We offer small and large pre-gelled Surface Electrodes designed specifically for EMG-recordings. They can be used as both active and reference electrodes. The large electrode can also be used as a ground electrode. The small electrode is especially suited for pediatrics.

- The electrodes are made with pre-gelled adhesives that are gentle on the skin yet keep the electrode firmly in place.

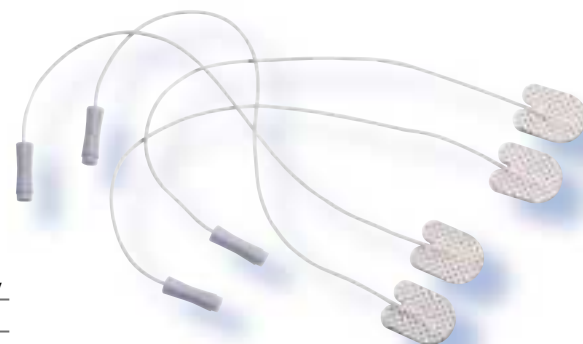


Item no.	Description	Size	Qty
9021S0231	Single-Use Pre-gelled Surface Electrode	Large	12
9013S0221	Single-Use Pre-gelled Surface Electrode	Small	12
9012C1161	EMG Cable for Surface Electrodes and Polygraf ID		1

### Sponge Electrodes

Our new Anal Sphincter and Vaginal Sponge electrodes have a patient-friendly design and offer several benefits:

- The sponge keeps the electrode firmly in place to eliminate movement artifacts.
- The slit in the sponge facilitates measurement of EMG and pressure simultaneously.

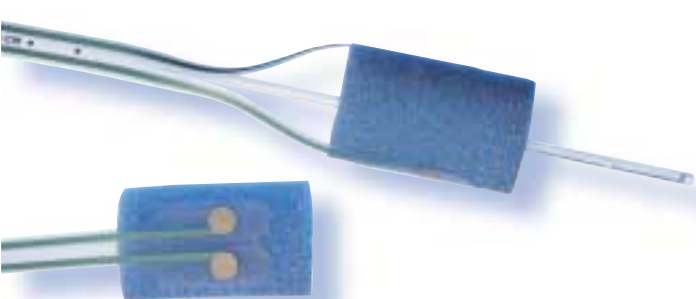

Item no.	Description	Qty
9013L4811	Single-Use Anal Sphincter Sponge Electrode	25
9021L0203	Single-Use Vaginal Sponge Electrode	25
9012C1451	EMG Cable for Sponge Electrodes and Polygraf ID	1

### Anal Plug Electrodes

The Anal Plug Electrodes allow you to measure the EMG signal of the anal sphincter muscle and come in two different sizes.

- The Anal Plug Electrodes feature a convenient pathway for inserting a rectal pressure catheter.



Item no.	Description	Size
9013K0782	Multi-Use Anal Plug Electrode	Large
9013K0792	Multi-Use Anal Plug Electrode	Small
9012C1201	EMG Cable for Plug Electrode and Polygraf ID	

### Ground Cables and Sensors

Item no.	Description	Qty
9021S0221	Grounding Strap Electrode	1
9013C0741	Ground Electrode Cable for Disposable Surface Electrodes	1